

The Nm1 (Newton Microscopes): Part 2 of 2

An in-depth examination and comparison to other folded-optics designs

R. Jordan Kreindler (USA)



Figure 1. Newton microscope in its foam lined plastic transport case as received

Introduction

The Newton Microscope is really three (perhaps one may say five or six) microscopes: the Nm1 - 400, Nm1 - 400 XY, Nm1 - 600, Nm1 - 600 XY, NM1 - 1000, and the Nm1 - 1000 XY. The primary differences are the magnifications available, and the presence or absence of a mechanical stage. Fig. 2 provides a close-up view of the Nm1 - 1000 with and without mechanical stage. On the Newton Microscope website, the microscopes are sometimes shown as Nm¹ rather than simply Nm1. Nm1 seems the more popular, and is the one used here.

At the time of this publication, the Newton Microscope site had not yet had a chance to post photographs of the Nm1 - 400 and Nm1 - 600 XY in their sales illustrations, and instead used pictures of the Nm1 - 1000, as a stand-in. All Newtons come with an external USB power cable and a "Rugged foam lined carry case". A picture of the inside of this case holding a Newton Nm1 - 1000 XY and accessories is shown in Fig. 1.

The Nm1 - 400, contains a 10x eyepiece and 2 objectives, 10x and 40x. The Nm1 - 600 XY has an additional 60x objective, and the Nm1 - 1000 XY has a 100x oil immersion objective in place of the 60x objective. If XY stages are present, these models contain an XY suffix after their designation. Thus, we have the Nm1 - 400 XY, the Nm1 - 600 XY, and the Nm1 - 1000 XY.

The Newton Microscopes prices are listed as Nm1 - 400 £399.95, Nm1 - 600 XY, £579.95, and the Nm1 - 1000 XY, £607.95. This may change over time, as more Nm1s are sold, or additional ways to reduce manufacturing costs are implemented. Newton prices have already changed since the author's purchase earlier this year. The mechanical stage can be purchased with the microscope as, for example, the Nm1 - 1000 XY in Fig 1 and Fig. 2b, or it can be purchased separately as the "XY Slide Indexer". When this article was published, the Indexer was listed at £79.95. In keeping with modern practice, most Newton microscopes and accessories have prices ending with £.95.

For some EEC buyers, and some others, final prices may include an additional tax on the purchase price, i.e., a VAT (Value Added Tax). For example, for UK buyers this tax is currently 20%, a rather high addition to the Newton Microscope prices quoted above. For other buyers outside England there may be other additional charges, e.g., customs duty.

Newton Microscopes-Nm1s



Figure 2. (a) Nm1 - 1000 with stage clips, and (b) the Nm1 - 1000 XY. Figure (a) Courtesy, and with permission, of Rick Dickinson

The mechanical stage, if obtained separately, comes in a plastic case, essentially identical in size to that provided to house the Meade Readiview, Fig. 3, and appears to be made of the same material.



Figure 3. Nm1 mechanical stage ("XY Slide Indexer") case

Newton Microscopes-Nm1s

Also available are adapters to connect various cameras to the Newtons. These include adapters for M42 mounts, £39.95, the Apple iPhone 4 and 5, £17.95, and Android phones, £17.95. Fig. 4 shows two Newton Nm1 camera adapters. Fig. 5 shows one of these adapters connecting the Nm1 and Olympus DSLR. In addition to the camera adapter, other accessories including a Cullman Nanomax 2000 tripod £34.95, and a Pelican Case £149 are available from the maker.

According to the Millennium-Microscope and Newton Microscopes sites, Newton Microscopes were initially developed to diagnose tropical diseases, e.g., malaria, tuberculosis, and helminthic diseases, i.e., diseases relating to parasitic worms, in third world countries, where relatively inexpensive microscopes were felt to be needed.

Owing to "real world" constraints the price of Millennium Health Microscopes (MHMs), [after early evaluation trials the name was changed to Newton Microscopes] were released at considerably higher prices than originally envisioned. To quote from a 2007 paper by Keith Dunning and J. Russell Stothard (Dunning, 2007), "*The MHM could be available within the next two years at no more than 75 USD.*" At today's exchange rate, an Nm1 - 400 costs approximately USD \$652.63, or almost 900% higher than the originally estimated price.

Considering that the Nm1 - 400 is now the only Nm1 version available for purchase without a mechanical stage, its price is relatively low. However, Nm1 - 400 users can purchase the optional XY Slide Indexer for £79.95. It is installed into the 400's slide clip stage openings. This makes the price of an Nm1 - 400 XY about £479.90, or almost £100 less than the cost of an Nm1 - 600 XY.

Newton Microscopes-Nm1s



Figure 4. Two Newton Microscope (Nm1) camera adapters

Newton Microscopes-Nm1s



Figure 5. Olympus DSLR mounted to Newton Microscope

Table 1, reproduced and updated from Part 1 of this paper, provides more details about the Nm1s.

<u>Newton Microscopes</u>					
<u>Name</u>	<u>Launch Date</u>	<u>Sales Date</u>	<u>~ Number Sold</u>	<u>Manufacturing Loc.</u>	<u>Launch Price</u>
Nm1 - 400	2013	2013	n/a	China	£412.50
Nm1 - 400 XY	2013	2013	n/a	China	£492.46
Nm1 - 600	2013	2013	n/a	China	£497
Nm1 - 600 XY	2013	2013	n/a	China	£579.96
Nm1 - 1000 XY	2011	2013	n/a	China	£607.46
Newton Fluorescence	unk	unk	n/a	China	est. £7-800

Table 1. Newton Microscopes

Overly optimistic forecasts for devices to be used in third world countries are, apparently, not uncommon. MIT's Dr. Negropte of the "One Laptop Per Child" project, established to develop \$100 laptops for developing countries, has had similar issues. That project's goals have, to date, proven elusive, as has the USD \$75 MHM. Low prices may often be based upon overly optimistic assumptions. These include, (1) forecasts that third-world governments, or other organizations, will buy a very significant number of units, thus allowing for possible cost savings through large-scale production, (2) prices for standard components used in these products will drop significantly, and/or (3) a failure to consider the other devices available to accomplish similar goals, etc.

For example, Newton Microscopes prices are comparable to many traditional binocular benchtop instruments, and these have expanded capabilities and options. Today, weight is often not the factor it once was, before the advent of reasonable roads and modern vehicles.

In spite of their higher than originally forecast prices, the Newton Microscopes and indeed all the Dickinson/Dunning folded-optics microscopes are unique. They are modern implementations of Dr. McArthur's folded-optics designs to three-dimensions. Newton Microscopes designed for third world applications deserve praise, in spite of being more costly than their original price goals, as it is hard to know what the prices for these instruments would have been had they had not originally targeted very low price points.

Design

The design goals for the Newton Microscopes were admirable. To quote from the "Introduction" section of their URL (newtonmicroscopes.com, 2013),

In contrast with infectious diseases in First World countries from viral and bacterial agents, diseases in the tropics largely result from a very different group of agents: the protists (single celled animals) and helminths (worms). ... entrenching many populations in grinding poverty... the greatest culprit is malaria which is present in 140 countries placing some 3.2 billion people at risk, while for helminths 650 million people are at risk of schistosomiasis.

Newton Microscopes-Nm1s

The Nm1s were designed to help identify these infectious agents, and thus reduce their impact. However, the Newton Microscope site also notes that while this was the original objective, these microscopes are also suited for use in a variety of other areas including general microscopy.

Newton Microscopes-Nm1s



Figure 6. Newton Microscope optical path. Courtesy, and with permission, of Rick Dickinson

Folded-Optics

All the Newton Microscopes use a unique compact 3D folded-optics design, see Fig. 6. Their size can be compared to other folded-optics microscopes available for sale on the used market, e.g., the McArthur, Vickers McArthur, the OU McArthur, the Nikon H, and the Swift FM-31. Also, new versions of FM-31 clones are now available for sale, (Kreindler, 2011-1). Comparison pictures, illustrating the relative sizes of the Newton Microscope in comparison to some other folded-optics models can be seen in Fig. 7. There were a variety of manufacturers who made metal McArthur microscopes. The two examples included here are typical. Other Dickinson/Dunning folded-optics microscopes were discussed in Part 1 of this paper. However, these models were not designed, as the Nm1 series was, or the other folded-optics microscopes shown (the OU McArthur being the exception) for serious professional use.



Figure 7. Newton Microscope and other folded-optics instruments. First Row: Nm1 - 1000 XY; Second Row: McArthur, Vickers McArthur, OU McArthur; Third Row: Swift FM-31, New FM-31 Clone, TWX-1, Nikon H

Newton Microscopes-Nm1s

The comparison Table 2, provides more details for the instruments shown in Fig. 7.

Comparison of Some Folded-Optics Microscopes⁽¹⁾			
Name	Approximate Size (mms)⁽²⁾	Approximate Weight (lbs,ozs)/grams	Comments
Newton Nm1 – 1000 XY	158 x 163 x 66mm	1lb, 6.2oz/629.4g	w mechanical stage
McArthur	101 x 50 x 88mm	1lb, 5.0oz/595.3g	w/o mechanical stage
Vickers McArthur	100 x 50 x 65mm	1lb, 2.0oz/510.3g	w/o mechanical stage
OU McArthur	129 x 25 x 75mm	----, 6.8oz/192.8g	w/o mechanical stage
Nikon H	150 x 65 x 109mm	1lb, 9.2oz/714.4g	w/o mechanical stage
TWX-1	131 x 40 x 115mm	1lb, 9.0oz/708.7g	w/o mechanical stage
Swift FM-31	147 x 115 x 108mm	2lb, 0.2oz/912.9g	w mechanical stage
FM-31 Clone	139 x 45 x 100mm	1lb, 12oz/793.8g	w/o mechanical stage

(1) All folded-optics microscopes were measured without their carry cases.
 (2) Length (front to back), width, height

Table 2. Comparison of some folded-optics microscopes

We can compare the Nm1 - 1000 XY's size and weight to the other folded-optics instruments presented here. It is larger in length and width, but shorter than almost all others in height, even with its mechanical stage mounted. Its mechanical stage is heavy, but four instruments in the table are heavier, including the Swift FM-31, also with a mechanical stage. The Swift FM-31 is almost 70% heavier than the Nm1 - 1000 XY, and for field portable microscopes weight is an important factor.

Unpacking and Use

My microscope was sent directly from the manufacturer. It arrived in a securely-taped corrugated cardboard box approximately 245mm x 100mm x 275mm. The box contained a foam lined plastic case with the Nm1 - 1000 XY. Inside that case was a small plastic case containing four slides: (Tilia X.S. of 3-year stem, Pollen Types W.M., Dicotyledon T.S, Sunflower, Zea mays, com L.S. of kernel [along with rice and wheat this is one of the top cereal crops worldwide]). The travel case also contained two Silica Gel packets, a fine focus adapter, an instruction manual, a USB cable, and stage clips, see Fig. 1. An adapter to connect the microscope to a camera was also shipped in the box, but was packed outside the plastic case.

Fig. 8 shows the foam lined plastic case closed, with microscope and manual on top, and included accessories in front. The fine focus adapter, an approximately 55mm diameter plastic component with a circular cutout that goes over the focus wheel, is not visible in the photo. It is shown dismounted and mounted in Fig. 9.

Out of the box, the microscope is almost ready to use. All that is required is that the fine focus adapter, if it is to be used, is mounted over the focusing knob. Reading the manual is helpful, but is not necessary to use the microscope, as most functions are easy to understand without instruction.

Newton Microscopes-Nm1s



Figure 8. The Newton Microscope, and accessories, outside of its shipping box. Slide on microscope is not included.



Figure 9. Newton Microscope fine focus adapter dismounted and mounted

Illumination

A slide is mounted on the stage, Fig. 8 and the illuminator is turned on via the button shown in Fig. 10. The power to the internal illuminator is provided by 3 AAA batteries located, and accessible, in a compartment underneath the microscope, Fig. 11. Illumination intensity is adjusted by a knurled knob, accessible from the underside or back of the microscope; it is shown in Fig. 12. The intensity control is graphically coded to show the direction of increased and decreased brightness.

Newton Microscopes-Nm1s

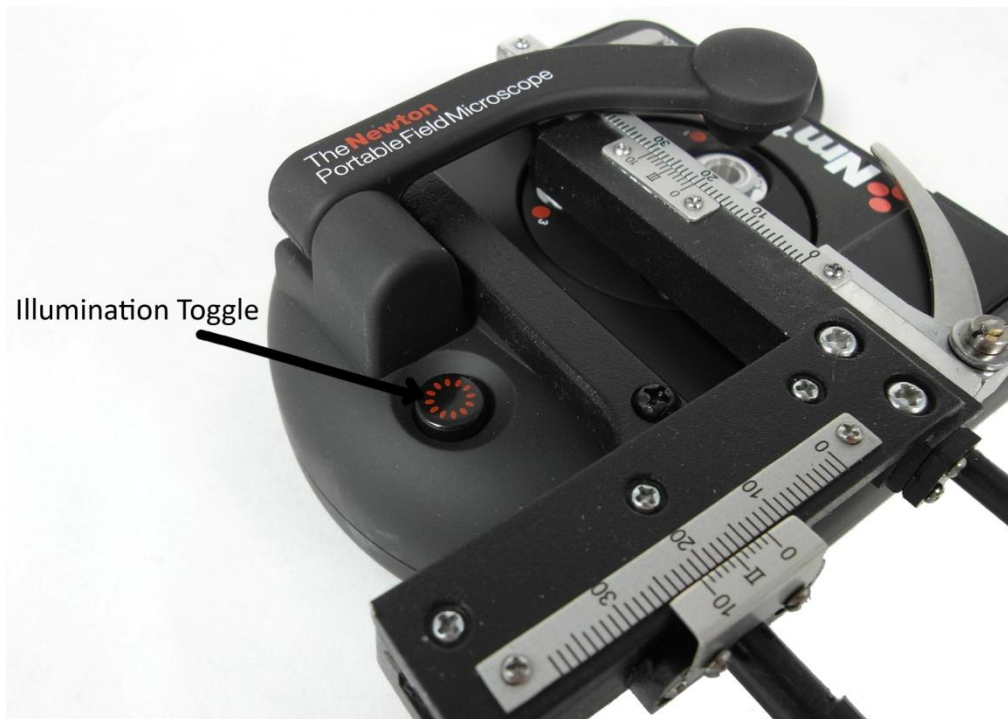


Figure 10. Newton Microscope illumination toggle

Newton Microscopes-Nm1s



Figure 11. Nm1 battery compartment with its three AAA batteries, underneath the microscope

Newton Microscopes-Nm1s



Figure 12. Knurled illumination adjustment knob

The built-in illumination arm (light arm) can be swung away, so that a slide can be lit by external illumination, Fig. 13. I found this most helpful, as the built-in white LED illumination was usually too bright for me at 10x magnification, even at its lowest brightness setting.



Figure 13. The light arm is lowered for built in illumination, and raised for external illumination

Eyepiece and objectives

The Nm1 -400, Nm1 - 600, and Nm 1- 1000 are all supplied with RMS 10x wide-field eyepieces. Newtons can have up to three objectives. The 400 comes with two, 10 x and 40x. The 600 with three, 10x, 40x, and 60x, and the 1000 uses a 100x oil immersion objective in place of the 60x in the Nm1 - 600. Thus magnifications run from 100x to, if the Nm1 - 1000 is selected, 1,000x. If the Nm1 - 400 is selected a blank plug is used in place of the third objective. Thus, all Newtons appear to be of the same design, except for the choice of objectives, and optional mechanical stage.

The field number (FN) of the Nm1s 10x eyepieces is about 18. This is lower than most benchtop instruments, which often have FNs of 20 and greater. However, the Nm1's eyepieces are certainly adequate for their intended screening purposes.

At its greatest, the human eye has an approximately 7mm opening. However, this only occurs when we are younger, and in very dark environments. As we age this opening contracts at the rate of about 1mm/year after our first three decades.

Newton Microscopes-Nm1s

The distance above the eyepiece where we can hold our eyes and see the full image field is referred to as eye relief.

On the Newton microscope site, the Nm1 is shown used in quite bright African sunlight. This means that the image from the microscope may be greater than the opening of the eye in these environments, i.e., the full optical capability of the microscope may not be applicable under those conditions. Of course, this would be true for any microscope used in the same circumstances, but less so for a microscope used on an inside benchtop, rather than outside in bright sunlight.

The Nm1 eyepiece was not designed to be used with eyeglasses, as its eye relief is too shallow. I have more expensive "High Eye-point" eyepieces on almost all my benchtop instruments, but these "High Eye-point" components are not available from Newton Microscopes. Users with astigmatism, over 50% of the population and over 80% of eyeglass users, would likely benefit from being able to keep their eyeglasses on when using an Nm1, so what they saw would not be, unnecessarily, distorted.

The objectives can be seen from the top of the microscope. Fig. 14 shows a view of the three objectives for the Nm1 - 1000 XY. The objectives, although they appear accessible, are not user exchangeable. Thus, e.g., if one desired a Newton microscope with 40x, 60x, and 100x, this would probably need to be specially configured by the manufacturer.

The turret and surround are black, which might result in heat issues in sunny climates.

Newton Microscopes-Nm1s

The small objectives are marked on their top surface, and their positions are marked on the turret. They use standard color codes (these codes essentially use the spectrum starting at black for the lowest magnification and ending with white as the highest, using cobalt blue, rather than violet, to identify 60 and 63x).



Figure 14. Nm1 -1000 objectives from top of microscope

Newton Microscopes-Nm1s

The color codes seen in Fig. 14 are 10x yellow, 40x light blue, and 100x white. The 60x objective is coded cobalt (i.e., a dark) blue. The immersion lens does not have a second color code to indicate its immersion media. The position locations on the Nm1 - 1000 are: 1 = 10x, 2 = 40x, and 3 = 100x. On the Nm1 - 600 the third position would contain the 60x rather than a 100x objective. The Newton Microscope URL, http://www.newtonmicroscopes.com/nm1_system.html, lists the lenses numerical apertures as 10x 0.25, 40x 0.65, 60x 0.80, and 100x as 1.25. This URL also provides detailed information about these microscopes.

A red dot outside the turret indicates the objective currently in use. In Fig. 14, it is the 40x.

The objectives are changed by a continually rotating knob, with click stops, on the bottom of the microscope, Fig. 15. This knob is graphically coded to show the direction of increased / decreased magnification.

Figure 15. Rotating knob under the microscope, used to select objectives



Newton Microscopes-Nm1s

The dimples (pin holes) on the objectives might suggest that they were there to allow users to change the objectives, e.g., to purchase the 60x and 100x and interchange these as needed.

However, to quote from an email to me from Rick Dickinson,

The design intent was precisely that - for the user to be able to add or remove an objective with a tool. This would not be a common thing and would probably be a swap between a 60 and 100, but still useful. However, the small scale of the objectives and the physical reality is such that it is extremely difficult, risky, and beyond the capability of most, and in some cases even rendered impossible due to mechanical effects. For example, I have a specially made tool to retrospectively insert a 60 or 100 into the blank position of a Newton 400, insertion is fine but extraction can be impossible if the objective is tightened down. Extraction of a factory fitted lens is impossible from outside due to the torque required to unlock the thread. Furthermore, the downward pressure required exerts mechanical leverage forces on the focus mechanism, which is undesirable. Then there is the concern for debris falling on the first mirror if the lens extract/insert is done over this location. Future production will remove the two 'pin holes' in the tops of the objectives.

—(Dickinson, 2013)

That is the objectives are not user changeable, and the current pin holes on the objectives will be removed in future releases of the Nm1.

Use

I found the Newton's unique almost circular shape, rather than the more traditional rectangular design, of a field portable microscope made it more difficult, for me, to hold comfortably. This impression may change over time. However, if the slide looked at will allow, the Newton can be turned on its side to look through the eyepiece, and this may prove a steadier way to hold the microscope for some users in the field. The eyepiece for the other field portable, folded-optics, microscopes shown in Fig. 7 all have eyepieces that look down into the microscope. The FM-31 and its clones are also tilted. I found these eyepiece angles, particularly those of the FM-31 and its clones, more comfortable to use. I also find it easier to use a field microscope with a flat bottom, with the eyepiece on top, as this allows me to place the microscope on a convenient surface for more stability when viewing.

However, the Nm1 series provides two tripod sockets beneath the microscope, located in both the front and back. The rear socket seems the appropriate connection if the microscope is to be used without an attached camera. A study tripod is required to produce stability, e.g., the Cullman Nanomax 200 T sold by Newton Microscope. With its compact ball head it can provide additional stability for the Newton. The 'T' here stands for travel, and this tripod weighs only about 700g. However, needing a tripod negatively affects portability. That is, it means extra bulk, and in these situations the Nm1 should be evaluated against more traditional instruments. Fortunately, the Cullman is relatively short even when fully extended, as it does not provide center to leg bracing. The Newton Microscope price for this tripod, £34.95, is less than the typical US price of almost \$60, and it can be recommended if you decide to purchase one of the Nm1s.

I found the microscope's tripod sockets somewhat shallow. I had a problem fully inserting the screw from a SLIK tripod with attached head, and could not lock a quick release plate from a Velbon tripod to the microscope. The knob to change objectives is just barely accessible, with difficulty, if a normal size tripod quick release is used, but presents no accessibility problem for a smaller size one. The shallow sockets should be taken into account if considering a tripod other than the one sold by Newton Microscopes.

For using the microscope when standing, a taller extending tripod with braces from the center column to the legs is probably desirable and many models are available, although not through Newton Microscopes.

If a camera is attached, the front socket may provide more balance. However, if a camera is used for documentation since the microscope is monocular, obviously the image viewed is not the same as that seen through the microscope, as it will be modified by the camera's optical system.

Optical Performance

Stability and usability are quite important but, possibly, the most important question for potential users is how the Nm1s optical systems perform compared to other folded-optics microscopes. For this analysis, I used the microscopes shown in Fig. 7 and Table 2 (excluding the OU McArthur). This table summarizes my subjective comparison of these microscopes in three areas: contrast, color fidelity, and sharpness. I did not attempt to show these factors in side-by-side photos, as I was more interested in how the microscope image would appear to a user in the field.

I have summarized my results in Table 3, where I used the following subjective rating scale for the three areas mentioned: Superb, Excellent, Very Good, Good, and Fair. I selected five slides covering a range of areas and objects for this analysis: two from the set supplied with the Newton and three others: (1) Tilia stem, (2) Pollen types, (3) Corpus ventriculi, (4) Jasper, Thin Section, and (5) Diatom mix.

I left the OU McArthur out of this summary, as it was not designed for professional field use. I was thus able to limit the resultant ratings, in Table 3, from Very Good through Superb.

Optical Comparison of Some Field Portable, Folded-Optics Microscopes

<i>Name</i>	<i>Contrast</i>	<i>Color Fidelity</i>	<i>Sharpness</i>
Newton Nm1 – 1000 XY	Excellent	Superb	Excellent-Very Good (1)
McArthur	Very Good	Very Good	Very Good
Vickers McArthur	Excellent	Superb	Excellent
Nikon H	Superb	Superb	Superb
TWX-1	Superb	Superb	Superb
Swift FM-31	Excellent	Superb	Excellent
FM-31 Clone	Excellent	Superb	Excellent

(1) It was difficult to get a consistently sharp focus using the Nm1's focus wheel even with the fine focus adapter. This may be this user's issue, or perhaps the focus gearing needs modification.

Table 3. Optical comparison of folded-optics microscopes

This quote, received from Mr. Dickinson, provides further information concerning the Nm1s optical performance,

It is theoretically correct to say that the single LED design does not allow the full resolution of the objectives to be attained, as the optimum lighting angle for each objective is not controlled, but in practical terms, it is difficult to see any difference. In fact for the 100x objective, which has an NA of 1.25, it is extremely rare that a conventional microscope is used correctly to gain such performance. The conventional Abbe condensers mostly have an NA of 1.25, but unless the condenser is oiled to the slide as well as the objective being oiled, the maximum resolution cannot be attained. Some of the more expensive aplanatic/achromatic condensers have NAs of 1.4, which means that most users waste even more of the performance (and money).

– (Dickinson, 2013)

Conclusion

For general microscopy, the Nm1 - 1000 XY was easy to use. The mechanical stage was a particular pleasure. It was smooth and easy to adjust, and its indexing capability let me return to a location of interest without problem. I also liked the swing-in and swing-out built-in LED lighting capability. It is possibly the easiest to use illuminator I have seen on a field portable microscope. Although at 10x it may provide too much light even at its lowest setting. I appreciated the open accessibility to the Nm1's stage, when the built-in LED is swung away. The only other folded-optics models I tested with similar easy accessibility to the stage when a mechanical stage is mounted are the no longer made Swift FM-31, and its current clones. The slide placement locations on later McArthur's, the Nikon H, and the TWX-1 restrict stage access.

The bottom-line question for potential buyers is probably, should I consider an Nm1 microscope. Based on my experience, particularly with its easy to use lighting system, its mechanical stage with graduated indexing, and its convenient stage access, I would say yes.

The only caveat, for me, is the Nm1's focusing mechanism.

I tested the Nm1 - 1000 XY's use in general microscopy, but not against its original design goal of detecting diseases in the third world. Since this was the original motivation for its development, one would hope its performance in this application would be exceptional.

However, Newton Microscopes felt the results of the tests, which they paid for, to be conducted at the Liverpool School of Tropical Medicine, were proprietary. They could not release names or contact information of the testers, in response to my request. Therefore, I could not directly learn the findings of these testers myself, and ask for answers to questions I, and perhaps readers here, would have liked to know about what the testers found.

If the built-in illumination mechanism of the Nm1 - 1000 XY, and its indexing mechanical stage are not important for you, I would recommend an FM-31 clone, as these can be found at less than 1/3 the price of an Nm1 - 400. The original Swift FM-31 was used extensively in the field with excellent results. The original Swift Company went bankrupt, and the new Swift Company no longer makes the FM-31. Fortunately, FM-31 clones are available new. However, the cover of the clone we examined is approximately 90% iron (Kreindler, 2011-1), so care must be exercised in wet and damp environments. Artificial light options are available as add-ons for the FM-31 and its clones, although these are not as easy to use as the illumination available on the Nm1s. FM-31 clones have a more traditional rectangular body, and eyepieces that are looked through at a downward angle. The rectangular body and eyepiece configuration means that they are relatively easy to rest on any convenient surface. See our earlier article in Micscape (Kreindler 2011-1) for a more detailed discussion of this microscope.

Optically the two best folded-optics microscopes I tested, in my subjective judgment, were the Nikon H and TWX-1. Of these, the ability to use slides positioned with the cover slip facing upward, made the TWX-1 my favorite. Both the Nikon H and TWX-1 produced superb images on all measures considered. In fact, they rivaled some monocular microscopes I have that were designed for benchtop use. However, both, as well as all the others instruments compared in Table 3 are monocular rather than binocular. This, for me, reduces their functionality significantly compared to binocular instruments, as I find seeing with both eyes allows me to detect features that might otherwise be difficult to locate.

I still find the Swift FM-31 a pleasure to use, and although it does not quite produce, for me, the quality of the images from the Nikon H and TWX-1, I still use it in the field as I find its mechanical stage, and tilted eyepiece, a pleasure to use. It uses mirrors, not prisms, and in spite of its aluminum alloy composition, as Table 2 shows, it is relatively heavy, especially with its mechanical stage attached.

The only microscopes in the Table available new today are the FM-31 clone and the Nm1 - 1000 XY. The clone tested did surprisingly well, but did not have an XY indexing stage. For new field portable microscopes, both the Nm1s and the FM-31 clones and can be recommended.

I would have liked to compare the results I found for the Nm1 - 1000 XY with the soon to be available Goren field portable microscopes (Goren, 2013), designed by Dr. Yuval Goren of Tel Aviv University (currently Visiting Lecturer University College London). Unfortunately, a production version of his design was not available at the time of this paper's publication. I hope these microscopes will be available in 2014. However, these do not use folded-optics. Although some should be available as binocular models, they are all likely to be larger instruments than the ones examined here.

If one can find used versions of the highly rated folded-optics instruments mentioned above (Nikon H, or TWX-1) in good working order, these should be considered as viable options, as should the original Swift FM-31.

It is nice to, once again, have so many new professional-grade models of folded-optics, field portable microscopes available.

Combined References and End Notes

Beale, M.A. (1991) *The Lensman portable microscope—an aid to school fieldwork?* Journal of Biological Education (Society of Biology) Summer 1991, Vol 25, Issue 2

Bracegirdle, Brian (1996) *Notes on Modern Microscope Manufacturers*. Oxford: Quekett Microscopical Club

Coghill, S.B. *Cytopathology*. 1992; 3(5):317-20

Dickinson, Rick (2013) Private emails. The author would like to thank Mr. Keith Dunning and Mr. Rick Dickinson for their valuable assistance in providing information and pictures. However, the content of this paper is the sole responsibility of the author

Dingley, Mike (1998) *The Lensman Microscope: a review*. Micscape Magazine
<http://www.microscopy-uk.org.uk/mag/art98/lensman.html>

Dunning, Keith (See Dickinson, 2013)

Dunning, Keith and J. Russell Stothard (2007). *From the McArthur to the Millennium Health Microscope (MHM): future developments in microscope miniaturisation for international health*, InFocus June 6, 2007

Enhelion (2002) *Micron by Enhelion*. Single sheet of instructions packaged with the microscope. July 2002

Flickr (undated-1)
<http://www.flickr.com/photos/9574086@N02/8757351327/sizes/l/in/set-72157624402374235/>, and
<http://www.flickr.com/photos/9574086@N02/8866715919/in/set-72157624402374235/>

Goren, Yuval (2013) *New Design for a Versatile Field Microscope*. Microscope Today, Volume 21, Number 3, May 2013

Journal of Microscopy (1996) Published on behalf of the Royal Microscopical Society.

Volume 183, Issue 2, pps 181-186, *Obituary for John Norrie McArthur 1901-1996*,

Wiley-Blackwell, August 1996. Available on-line at

<http://onlinelibrary.wiley.com/doi/10.1046/j.1365-2818.1996.1050648.x/abstract>

Kreindler, R. Jordan and Yuval Goren (2011-1) *Comparison of the Swift FM-31 Portable Field Microscope and an FM-31 Clone*. Micscape Magazine. March 2011

http://www.microscopy-uk.org.uk/mag/artmar11/A_Comparison_of_FM-31_and_clone.pdf

Kreindler, R. Jordan and Yuval Goren (2011-2) *The TWX-1 Folded-Optics Microscope: Monarch of the Folded-Optics Kingdom*. Micscape Magazine. November 2011

<http://www.microscopy-uk.org.uk/mag/artnov11/TWX-1.pdf>

newtonmicroscopes.com (2013). The Newton Microscope site, describing the Newtons and their applications, with a section to purchase these instruments

Nikon (2013) Model H Microscope, MicroscopyU, The Source of Microscopy Education.

<http://www.microscopyu.com/museum/modelh.html>

Rowntree, Frank (1996)

<http://www.microscopy-uk.org.uk/mag/articles/obituar1.html>

Saitoh, Motohiko (2006) McArthur microscope and Japan,

<http://www.geocities.jp/saitohmoto/essay/mcarthur/mcarthur-en.html>

The Herald, Glasgow (1996) *Obituary of Dr John Norrie McArthur, MRCS, LRCP, Hon FRMS, Hon Doctorate Open University, Hon FRS, born Glasgow, November 9, 1901, died Cambridge, April 26, 1996*. May 11, 1996

Walker, David (2005) *Microscopy on the move: A look at the Trekker field microscope*.

<http://www.microscopy-uk.org.uk/mag/artjun05/dwtrekker.html>

Walker, Spike (2013) Photomicrographer Spike Walker was kind enough to share his recollections of Dr. McArthur and his work with the Lensman organization, in private correspondence with the author. [Spike Walker is no relation to David Walker, mentioned just above.]

©2013 Text and photographs (except as noted) by the author.

The author welcomes any suggestions for corrections or improvements. He has a continuing interest in folded-optics and field portable microscopes, and stereo and unique historical microscopes.

The author can be reached at:

R. Jordan Kreindler: leona111rjk@gmail.com

Published in the online magazine Micscape, December 2013

www.microscopy-uk.org.uk

Please report any Web problems or offer general comments to the Micscape Editor.

Micscape is the on-line monthly magazine of the Microscopy UK web
site at
Microscopy-UK