

Newton Microscopes-Nm1s

Folded-Optics Microscopes 1 of 2

The Nm1 (Newton Microscopes): Their Heritage

R. Jordan Kreindler (USA)

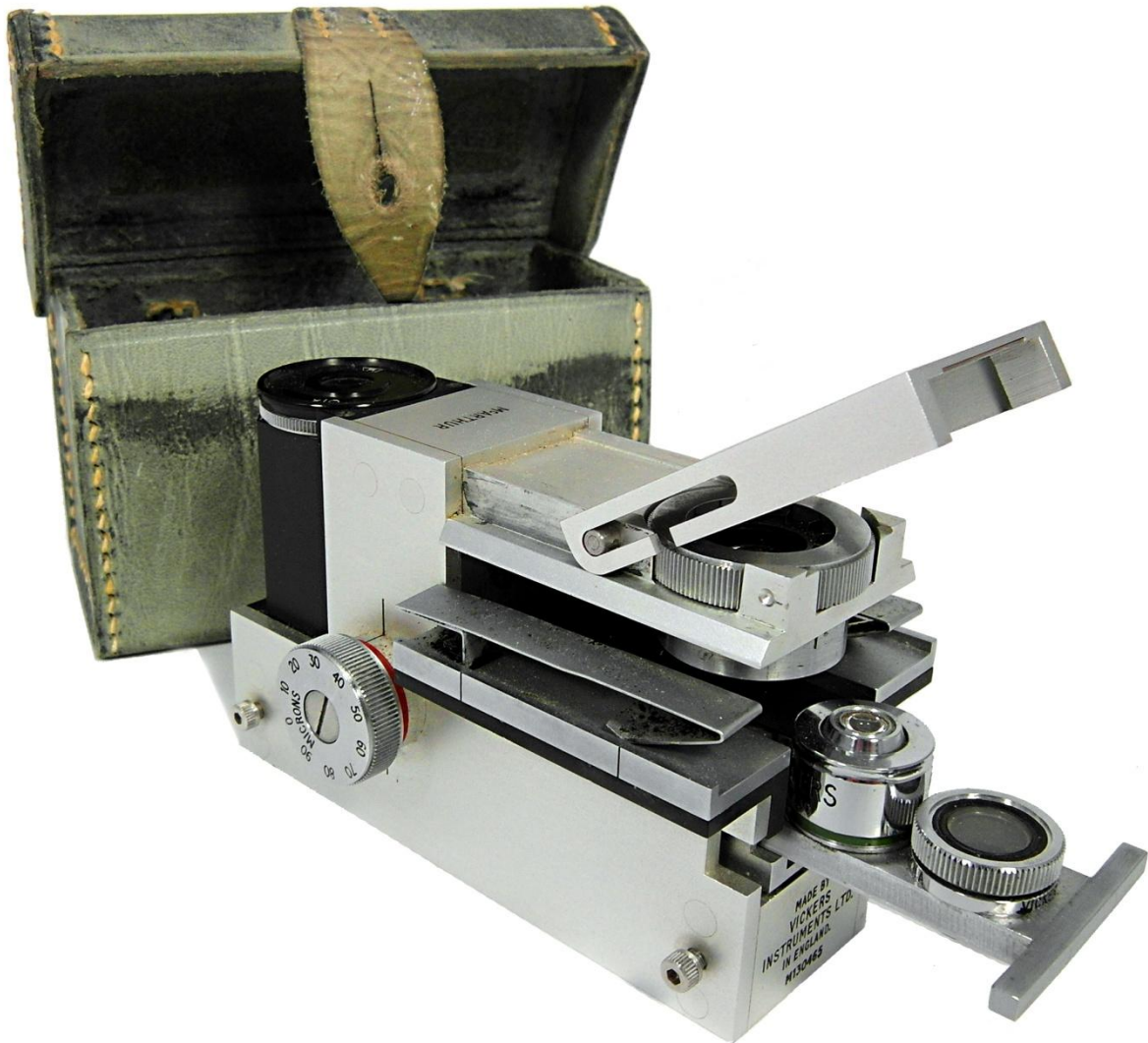


Figure 1. McArthur folded-optics microscope, made by Vickers Instruments, Ltd. England.
(In case, approximate size: 11cm x 7cm x 6cm; weight: 1 pound, 5.6 ounces). Post 1963

Folded Optics Microscopes

At the beginning of the 1930s a seminal, and for field microscopes a "sea change", event took place. Dr. John Norrie McArthur (1901-1996) invented and produced the prototype of his first folded-optics microscopes. He discussed his microscope before the Quekett Microscopical Club in 1933. McArthur's name is now synonymous with his folded-optics microscope. This seems fully appropriate considering his importance in the evolution of portable field microscopes.

His compact instrument was approximately 10 cm x 6.5 cm x 5.0 cm. During WW II, Dr. McArthur was captured by the Japanese in Borneo, but allowed to continue his research by Japan's chief civilian administrator for Borneo, Tatsuo Inagawa-san. This freedom allowed Dr. McArthur a chance to use his newly designed microscope. This experience suggested to Dr. McArthur that some important changes to his original folded-optics design were needed. He made these changes at the war's conclusion. (Kreindler, 2011-2)

Folded-optics microscopes can trace their lineage directly back to the seminal work of Dr. John Norrie McArthur, begun in the 1930s. According to the 1974 Guinness Book of Records, Dr. McArthur's microscope was the world's smallest high power microscope. (Saitoh, 2006).

Mr. Saitoh was a Japanese physicist who was fortunate enough to meet Dr. McArthur in 1972/73 and then a decade later in 1982/83. The following is from Mr. Saitoh's personal account, with minor English corrections by the author, as told to Mr. Saitoh by Dr. McArthur,

I was studying malaria in Borneo, and was captured and confined in a camp by the Japanese army during World War II. Fortunately, however, Mr. Inagawa, a high ranking officer of the army administration, understood the importance of the research, and encouraged me to continue. He provided me with supporting money for research and even milk for my baby, under the difficult circumstances of getting foods.

During my research, I desperately felt the need for a portable microscope usable in fields and under the kerosene lamp in a tent. After my release from imprisonment, I started developing compact microscopes with full versatility suitable for research. I am currently very busy in pursuing the ultimate design of models of various types, including the surface inspection model, to my complete satisfaction. (Saitoh, 2006).

Newton Microscopes-Nm1s



Figure 2. Plastic body McArthur microscope made for the Open University (OU), c. 1970s

The ubiquitous Open University (OU) McArthur microscope, Figs 2 and 3, perhaps best exemplifies McArthur's folded optics design. This microscope was made by Scientific Optics Ltd., Ponswood Industrial Estate, Hasting, Sussex England, for the OU; it was adapted for this use under Dr. McArthur's guidance. Traditional McArthurs (e.g., Figs. 1 and 4) are made of metal, and their cost was prohibitive for many users and uses. The unique OU McArthur

Newton Microscopes-Nm1s

microscope was made of plastic (acrylonitrile butadiene styrene) and was, relatively, inexpensive to manufacture.

However, using inexpensive construction materials has had its effects on long-term durability. Now, decades later, many of the OU McArthurs show a significant loss, if not complete lack, of functionality. However, in spite of the, arguably, substandard construction material, particularly in making the objectives slider, possibly the instrument's "Achilles heel", many OU McArthurs are still in use. This probably demonstrates the general robustness of McArthur's design, even when less desirable construction materials are used. Many OU models have been converted to use AAA batteries and LED illumination, as the original battery compartment was too narrow to accommodate AAAs without modification.

Note the warning in Fig. 2, to "SLACKEN Focus screw after use". This screw was used to focus the instrument by flexing the objectives carrier. If left tightened, the plastic can lose its flexibility and fixed in place would be non-useable for focusing. This is one reason many OU models are now non-functional, and why the objective slider is likely the "weakest link" in the microscope's design.

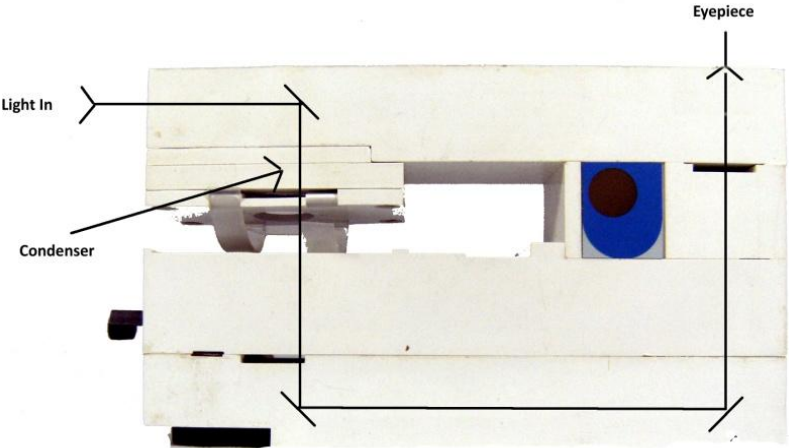


Figure 3. The folded-optical path in an OU McArthur microscope

As can be seen in Fig. 3, external light entering the McArthur is "folded" (reflected) through three right angles. It is first reflected downward to a second reflector, where it is transmitted across the base of the microscope to another reflector, and then reflected upward into the eyepiece. A slider can be moved so that the external light source is replaced by illumination from an internal light, e.g., lamp, or LED if modified.

Newton Microscopes-Nm1s

Objects are seen "upside down" compared to conventional microscopes. The object to be viewed faces downward, as in an inverted microscope, rather than upwards toward the eyepiece. Because of this "inverted" nature of the microscope, only a single fine focus is needed, as the object viewed will always rest on top of the platform, i.e., at a almost constant distance from the top of an objective.

In less expensive McArthurs, as in some recent SLR cameras, and many folded-optics descendants, e.g., the FM-31, mirrors are used in place of prisms to reduce manufacturing costs.

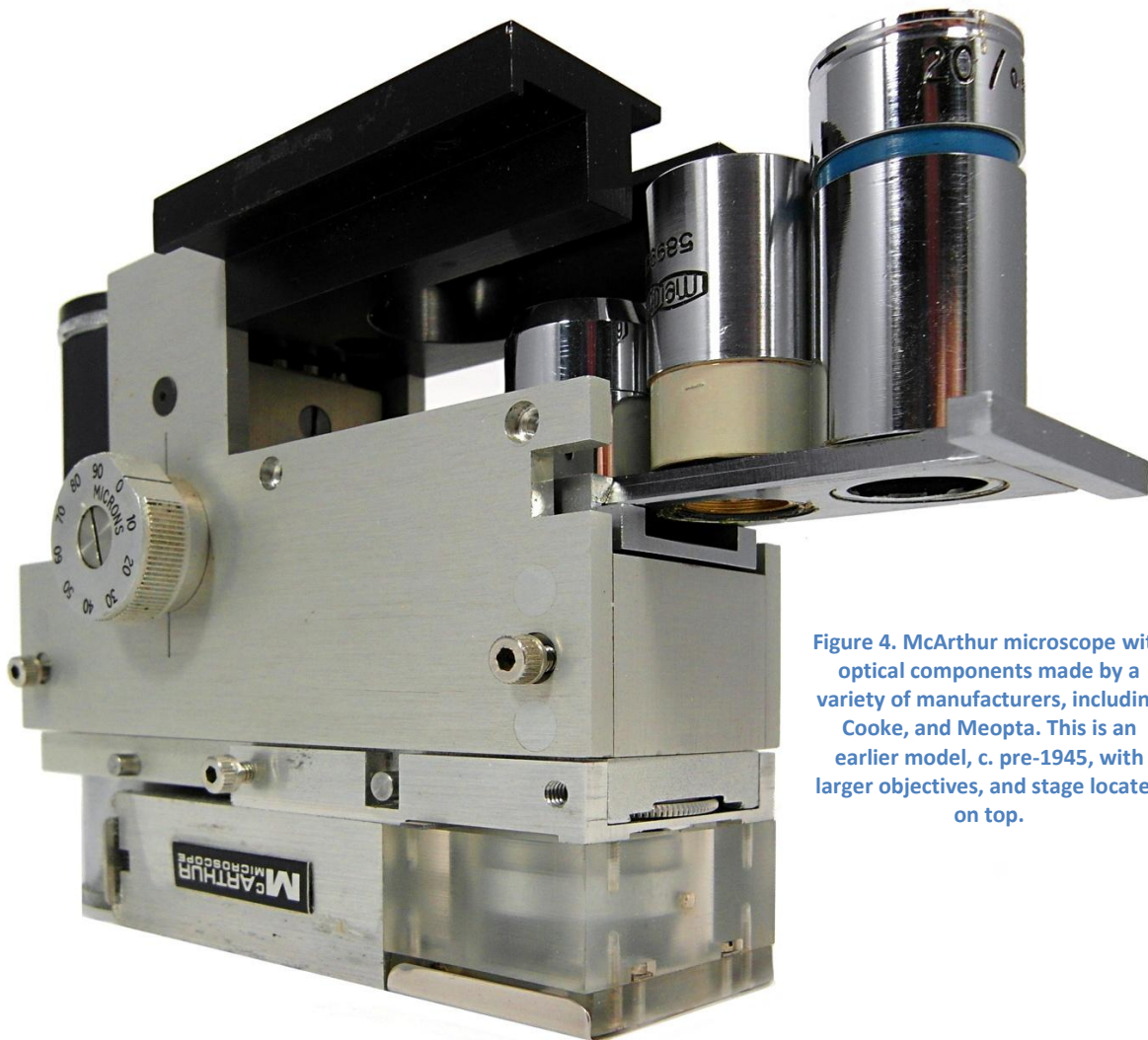


Figure 4. McArthur microscope with optical components made by a variety of manufacturers, including Cooke, and Meopta. This is an earlier model, c. pre-1945, with larger objectives, and stage located on top.

Newton Microscopes-Nm1s

McArthur's folded-optics design spawned a number of descendants. These include the Nikon H, Fig. 5, (Nikon, 2013),



Figure 5. Nikon H folded-optics microscope, c. 1960s

Newton Microscopes-Nm1s

and the TWX-1, Fig. 6, (Kreindler, 2011-2), as well as Tiyoda Models VI and VII, and.

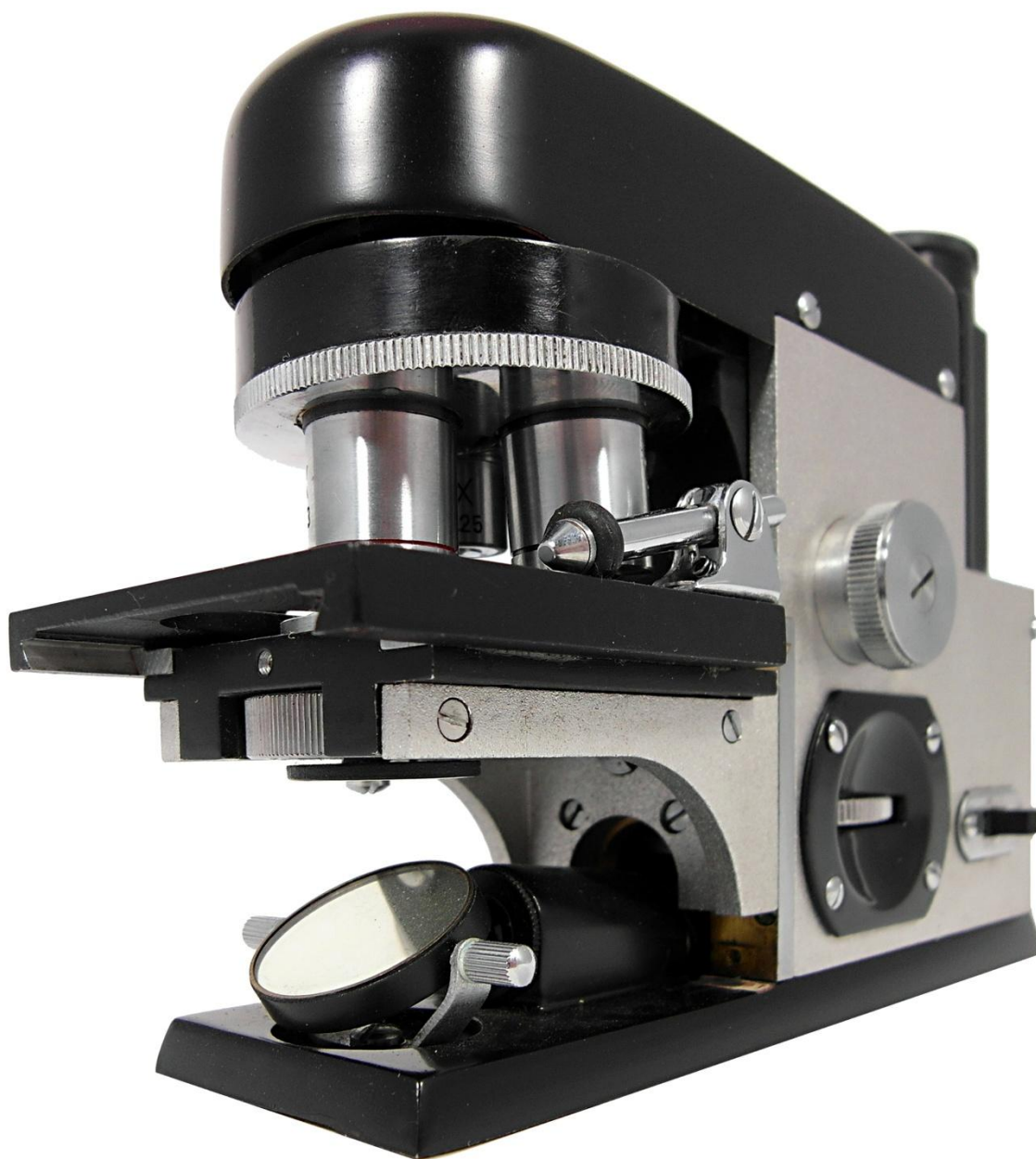


Figure 6. TWX-1, c. mid-1970s

Newton Microscopes-Nm1s

the Swift FM-31, Fig. 7,



Figure 7. Swift FM-31 with Maglite® illumination c. 1980s

Newton Microscopes-Nm1s

and its clones, Fig. 8, (Kreindler, 2011-1).



Figure 8. Current FM-31 clone with attachable illuminator, stored in original packing Styrofoam

Newton Microscopes-Nm1s

These descendents also include the folded-optics implementations of Keith Dunning and Rick Dickinson, such as the Lensman, Enhelion Micron, Meade ReadView, Trekker, and the Newton microscopes (the Newtons are discussed in detail in the second Part of this paper). To date, these Newton folded-optics microscopes include the Nm1-400, Nm1-400 XY, Nm1-600, Nm1-600 XY, and Nm1-1000 XY. Newton microscopes were originally marketed by Cambridge Optronics, but the name of the marketing/sales arm for Newton microscopes recently changed to "Newton microscopes" to avoid confusion. Cambridge Optronics will remain as a holding company, but "Newton microscopes" will be the public face of Cambridge Optronics for Newton microscopes.

The following table provides summary information on the folded-optics microscopes developed by Messrs. Dickinson and Dunning (Dickinson, 2013).

Table 1

| <u>Name</u> | <u>Launch Date</u> | <u>Sales Date</u> | <u>~ Number Sold</u> | <u>Manufacturing Loc.</u> | <u>Launch Price</u> |
|---------------------------|--------------------|-------------------|----------------------|---------------------------|---------------------|
| Lensman | 1988 | 1989 | 7,000 | UK | £99 |
| Enhelion Micron | Modified Lensman | 2001 | 4,500 | China | £99 |
| Meade ReadView | Modified Lensman | 2003 | 5,000 | China | £ 90 |
| Trekker | 2003 | 2004 | 7,400 | China | £ 49 |
| <i>Newton Microscopes</i> | | | | | |
| Nm1-400 | 2013 | 2013 | n/a | China | £ 412.50 |
| Nm1-400 XY | 2013 | 2013 | n/a | China | £ 492.46 |
| Nm1-600 | 2013 | 2013 | n/a | China | £ 497 |
| Nm1-600 XY | 2013 | 2013 | n/a | China | £ 579.96 |
| Nm1-1000 XY | 2011 | 2013 | n/a | China | £ 607.46 |
| Newton Fluorescence | unk | unk | n/a | China | est. £7-800 |

Newton Microscopes-Nm1s

Pictures and brief discussions of the Lensman, Enhelion Micron, Meade ReadView, and Trekker are provided below. These microscopes were designed as relatively inexpensive, highly portable instruments that could be carried into the field, and they met these goals quite well. They were not designed to meet demanding scientific requirements. Although they can be used for initial field or bedside screening, further definitive, detailed analyses requires a higher quality microscope with greater resolution. The Lensman and its derivatives are unique in using 3D folding, instead of the usual 2D folding found in most other folded-optics McArthur derivatives. They use a third reflective surface to fold (bend) the light in the z-dimension. These unique designs allow for a relatively large stage that is unencumbered by intrusions, and thus provide room for a removable ring magnet, visible in Figs. 11 and 12, to hold slides to a metal stage.

The Lensman, and its cohorts



The Lensman, Fig. 9, is roughly 4" in diameter; it was sold by Science of Cambridge, Cambridgeshire. It offers two magnifications, 80x and 200x. However, as with the Enhelion Micron discussed below, changing magnifications requires refocusing, as the dual magnification is achieved by moving a Barlow lens into the optical path.

Figure 9. The Lensman original concept. Illustration courtesy, and with permission, of Rick Dickinson

The optics for the Lensman were designed by UK optical designer Eddie Judd. In addition to dual magnification options, there are dual selectable illumination settings. The collapsible arm contains the bulb, which is powered by two AA batteries located in the saucer like section of the microscope, and accessible from beneath, Fig.9.

Newton Microscopes-Nm1s

Mike Dingley (Dingley, 1998), notes that the "*quality of the images are at best reasonable but not good enough to take photomicrographs ...*". M. A. Beale comments (Beale, 1991) that "*The Lensman is not a high-resolution piece of equipment, and consequently the images seen are poorer than with conventional microscopes. Focusing is rather difficult because the focus wheel is not locked tightly into place: when pressure is applied to the wheel, the image clicks out of view. ... At present the Lensman is an educational toy, and it has little value as a teaching aid. "*

In spite of these somewhat negative appraisals, the Lensman was unique and garnered a number of awards. These included the 1990 Prince of Wales Design Award, the 1990 BBC Award, the 1991 Archimedes Award, and the 1991 Camera of the Year Award, in the accessories category.

The Lensman could be purchased as a package with three slides, a 26 page booklet (the instruction manual was written by Professor Brian Ford), an instruction poster, and material on butterflies and moths, and other insects. This material was presented on microfiche, hence the square field of view. The use of microfiche (called micropage in Lensman "speak") was felt to make it easier for identification in the field without the need to carry field identification books.

However, this use of microfiche never caught on. Perhaps this was due to the rapid advances in electronic imaging that made microfiche seem "old-fashioned". Whatever the reason, the microfiche identification guides did not continue into subsequent microscope models. However, the more expensive square viewing format did persist, probably for its benefits to photography, and possibly the value of not having to retool. The numerical aperture (NA) for the Lensman and its derivatives is thought to range between 0.10 to 0.12.

A green colored ("Green Lensman") was made for the natural History Museum, London (Flickr, undated-1).

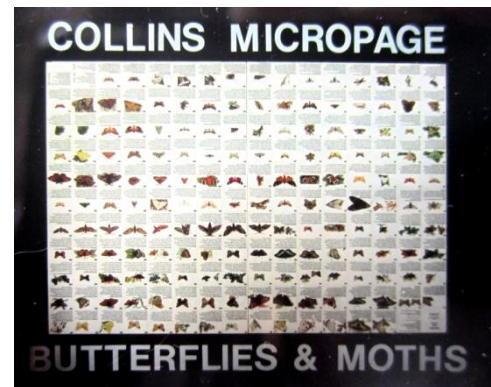


Figure 10. Lensman microfiche (micropage) used for field identification. Courtesy, and with permission of Rick Dickinson

The Enhelion Micron

The Enhelion Pro 160's one sheet guide in introducing the Micron states,

The Micron range of microscopes, designed and produced by Enhelion is a radical and innovative departure from microscope convention. They are designed with portability and ease of use as a priority offering a compact and versatile instrument solution with an extremely wide range of applications. ... (Enhelion, 2002)

Fig. 11 shows the Enhelion Micron in its storage case, in front of its retail box. This microscope provides both 80x and 160x magnifications, accessible via a knurled, rotating knob. However, the two magnifications still required re-focusing owing again to the use of a Barlow lens, i.e., they were not parfocal. The Micron uses two AA batteries for illumination. This is a light, highly portable instrument weighing, in its case, approximately 12.8 ounces. My microscope kit came with two slides: a somewhat amateurishly prepared slide of a "Butterfly Leg", and a clear shallow well slide. This microscope has a "Nextel" paint finish, and laser etched graphics.

The Meade ReadView

Meade ReadView's single user instruction sheet says the following,

ReadView (Fig. 12) is a radical new microscope design. Designed for portability and ease of use. ReadView is highly versatile and can fulfill a wide range of applications. ReadView provides two magnifications levels, 80x and 160x. ... The eyepiece may be removed (just gently pulled from the body of the microscope) and be used as a 10x viewer (Fig. 13).

ReadView has a lighting arm, which conveniently folds down and out of the way for storage. Power for the lighting arm is provided by two user-supplied AA-size batteries. Included with [the] kit is a magnetic holder to secure your slides or specimens in place, plus two slides: one prepared with a specimen and one left blank so you can prepare your own specimens.

Newton Microscopes-Nm1s

Enhelion was contacted by Meade to see if they would be interested in making the Enhelion Micron Pro 160 for Meade, with the "new" product marketed under Meade's name. Enhelion agreed. The primary change from the Enhelion model was the introduction of LED illumination, then under development for the Trekker (see below).

There were also some small cosmetic modifications such as the color of the microscope's outside diameter band and eyecup, the ReadView name radiating inward toward the origin of the microscope's circle, and the Meade name on the circumference. A newly designed retail box and new instructions were produced by Meade. The ReadView microscope was produced for Meade by Enhelion, in the same Chinese facility where the Micron was produced.

The Meade ReadView, as noted, is a slight modification of the Enhelion Micron, which in turn was a slight modification of the original Lensman. This microscope, as the one discussed above, comes with a prepared slide of a "Butterfly Leg", and a clear shallow well slide. Although, relatively inexpensive, the optics are glass not plastic. The ReadView's appearance is similar to that of the Enhelion Micron, Figs. 11 and 12.

Mr. Dickinson was kind enough to provide the following additional information (very slightly reformatted by the author) concerning the Lensman, Enhelion Micron, and ReadView.

- *All came from the same tooling.*
- *All had the same square field of view.*
- *All have an internal polished metal reflector to direct illumination upwards for solids viewing when the lighting arm is in the down position and illumination left on. (This aspect seems to be often missed, and was a useful function for hobbyists.)*
- *Slightly angling the lighting arm can introduce dark ground illumination effects.*
- *All microscopes had pull out eyepieces, allowing them to be used as a 10x loupe.*
- *To save cost on the optics, mechanisms, and space, all had a single objective lens and the dual magnification was achieved by swiveling a Barlow lens into the optic path. This is also why the two magnifications were not parfocal.*
- *All instruments were supplied with a molded magnetic ring to hold down specimens or slides on their ferrous (magnetic attractive) stages.*

--- (Dickinson, 2013)

Newton Microscopes-Nm1s



Figure 11. Enhelion Micron

Newton Microscopes-Nm1s



Figure 12. Meade Readiview handheld microscope



Figure 13. Readiview eyepiece removed for use as a 10x magnifier

The Trekker

The Trekker, Fig. 14, was purchased in a see-through case, with the eyetube of the microscope extended to the graphics on the box. It comes with an attractive drawstring, gray and yellow cloth pouch case, approximately 6-1/2" x 4". A hard-sided carrying (storage) case is also available. That case is similar to those provided for the Lensman derivatives discussed above. The author's new Trekker came with a prepared "Feather" slide.

One small, but meaningful, change was to the magnetic slide holder, which now had a matte aluminum-toned finish, making it harder to lose in the field.

The following comes from the back of the Trekker retail box,

- *High resolution coated glass optics*
- *Views solid and transparent objects and fluids*
- *Magnifies x35 and x10*
- *Unique square field of view*

The 10x specification, although accurate, may be somewhat misleading. This magnification is not obtained through the compound microscope, but by removing the eyepiece and using it as a separate magnifier. While handy, a separate 10x magnifier would be smaller, and provide more flexibility and usability than the removed eyepiece.

The approximate weight for the microscope with slide holder, retail box, and cloth case is 10.2 ounces. The Trekker with only the magnetic holder weighs approximately 6.2 ounces. The hard-sided case adds an additional 6 ounces.

In this foursome of folded-optics microscopes, the Trekker is somewhat the "odd man out". While it contains some parts from the Lensman family, it was/is an essentially newly produced variation of the Lensman. The re-tooling was funded out of an upfront royalty payment from Meade. Owing to the lower magnification, i.e., 35x, the lighting arm present on the earlier three models was no longer needed. Rather a daylight temperature, white light LED is placed inside the body to illuminate solid objects, when appropriate, from the side. In some instances, this LED placement can show more object detail, although occasionally at the expense of not adequately illuminating object areas with relatively deep cavities. The object viewing circle is surrounded by a transparent acrylic ring to diffuse the LED light more widely.

Newton Microscopes-Nm1s

The LED illumination required a different voltage than the earlier lamp's AA batteries, so provisions for a standard rectangular 9-volt battery were provided for in the battery compartment on the underside of the microscope, Fig. 15. Using a quality battery, illumination should be available for almost two weeks of continuous 8 hour/day use, without the need to change batteries.

In June 2005, David Walker evaluated the Trekker microscope for Micscape Magazine, after purchasing a new one for about 30 GBPs on eBay UK.

His article is interesting reading, as it provides both pictures, one with the microscope partially disassembled, and an independent evaluation of Trekker functionality. Mr. Walker concludes,

This is an elegantly designed microscope, genuinely pocketable and gives informative images of suitable subjects. It's also good value at its discounted price of ca. £30. The 35x mag would perhaps best suit a naturalist seeking a higher mag than a hand lens can offer to study the smaller features of fauna and flora critical to identification in the field ... It's also useful for studying aspects of smaller aquatic and land invertebrates. ... The Trekker is less suitable for more general studies of fauna and flora in the field where the lower power of a 5x - 10x hand lens with its greater depth of field, field of view and more intuitive use would be better. ... The Trekker's eyepiece does act as a 10x lens and is a useful feature, although ... its short operating distance makes it a less versatile hand lens than a dedicated one.

--- (Walker, 2005)

Mr. Walker also notes, *at 35x the depth of field is inevitably quite shallow (< 1mm)*, and further observes that the microscope is best with relatively flat subjects.

A T-mount SLR camera adapter, that attaches to the microscope on one side and accepts an SLR's T-mount adapter on the other, is available. The T-mount option assures compatibility with almost any SLR. However, photographic images will not be up to benchtop quality, as resolution is likely at best to be about 4 microns.

For me, of the Dickinson/Dunning microscopes discussed here, the Trekker, possibly because of its reduced 35x magnification, has proven the easiest to use.

Newton Microscopes-Nm1s



Figure 14. Trekker

Newton Microscopes-Nm1s



Figure 15. Trekker battery compartment (bottom) and object lighting (top)

Newton Microscopes-Nm1s

The Trekker is easily mounted on a tripod, Fig. 16, using the tripod socket on the underside of the microscope, shown in Fig 15. This allows for steadying the microscope for the use, or for photographs through the instrument.



Figure 16. The Trekker on a tripod with slide held by magnetic slide holder

Newton Microscopes-Nm1s

A picture of the "feather" slide provided with the Trekker is shown in Fig. 17. As can be seen, the resolution is not high, but likely acceptable for initial screening use in the field, and for hobbyists.

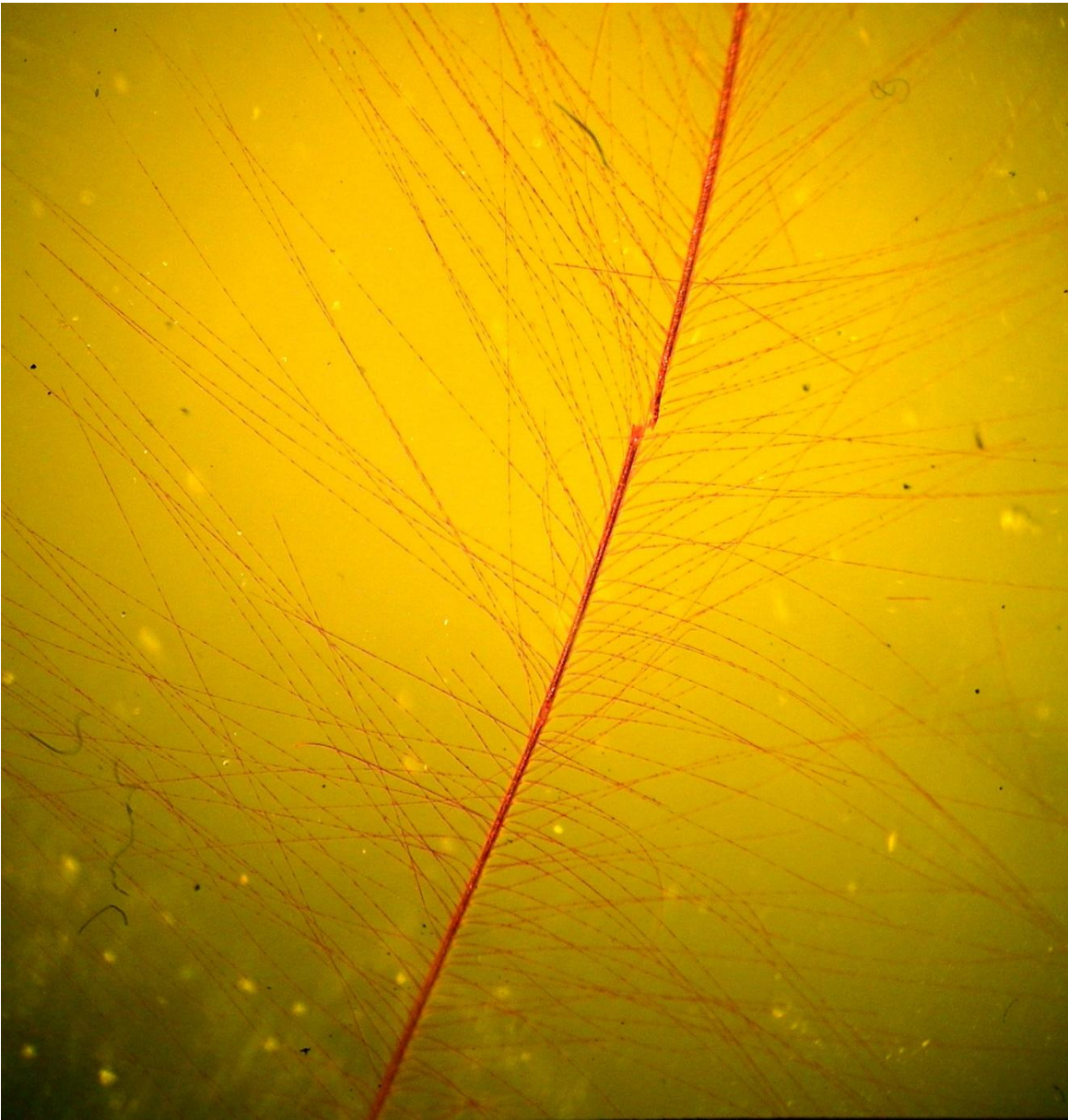


Figure 17. Feather slide provided with Trekker

Newton Microscopes-Nm1s

The Trekker evolved from the Enhelion Pico, Pico Peep 35, prototype platform, Fig. 18. Note the lack of a lighting arm that is present in the Pico prototype, but absent on the Trekker.

At roughly this time Messrs. Dickinson and Dunning started a new brand called "Looksmall". The website, www.looksmall.com, will be found printed on the underside of the Trekker.



Figure 18. The Trekker and its development predecessor the Pico. Courtesy, and with permission, of Rick Dickinson

Newton Microscopes-Nm1s

There were thoughts of a stereomicroscope in the Lensman family, but that, sadly in the author's opinion, never went beyond the concept stage, Fig. 19.

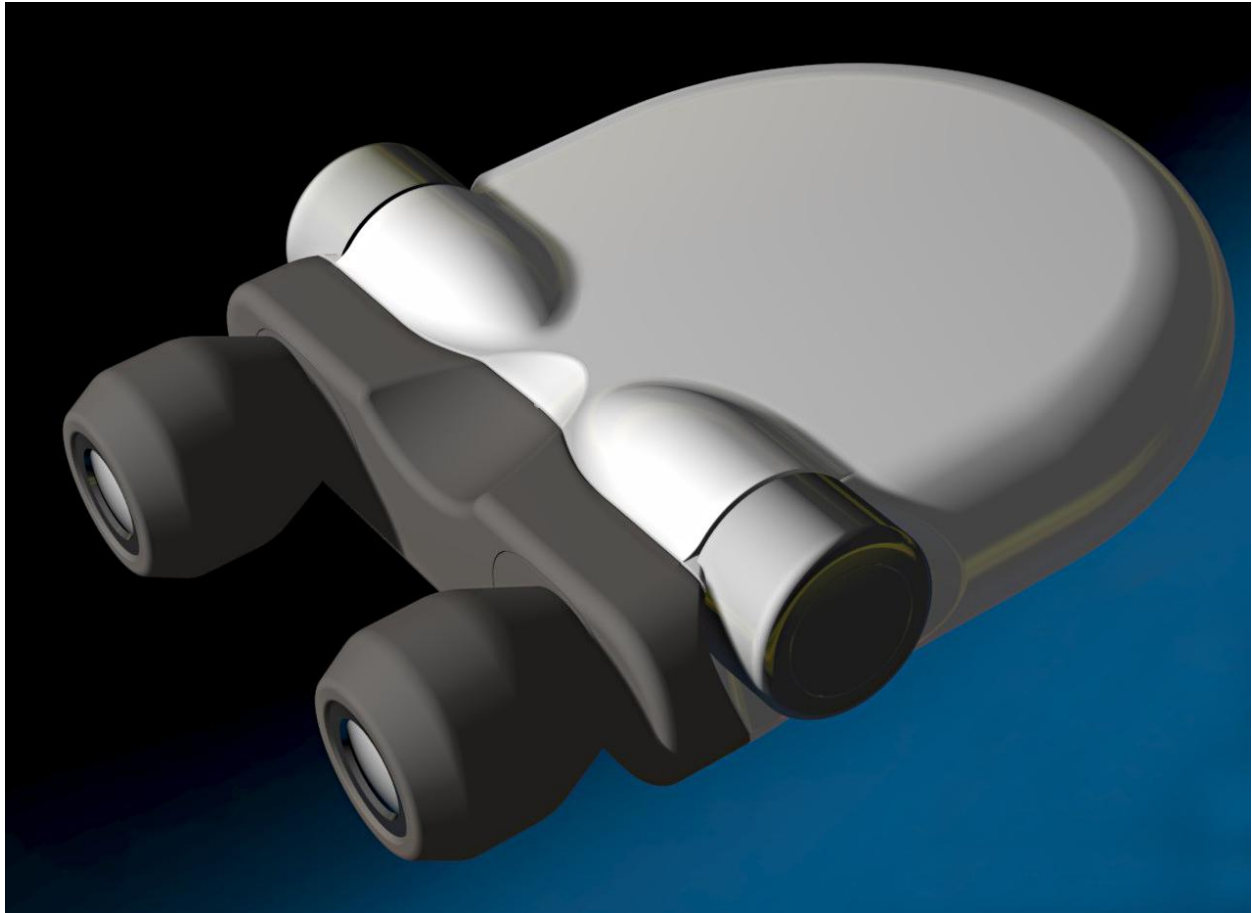


Figure 19. Concept design only for a stereomicroscope in the Lensman family. Courtesy, and with permission, of Rick Dickinson

Newton Microscopes-Nm1s

In the second Part of this paper, the Newton microscopes, Nm1-400, Nm1-400 XY, Nm1-600, Nm1-600 XY, and Nm1-1000 XY will be discussed in more detail. These microscopes are derivatives of the microscopes, discussed above.

These new Dickinson/Dunning folded-optics microscopes are now targeted toward a user community of scientists and other professionals. The Newtons have considerably increased capabilities compared to, e.g., the Lensman and its derivatives. Not surprisingly, a side effect of this increased capability is increased manufacturing costs. This has resulted in the Newton microscopes costing considerably more than their less capable predecessors, and exceeding their originally estimated prices.

As the Newtons are now targeted to a different audience than the Lensman and its descendants, it is appropriate to contrast them with, arguably, competitive professional-level field portable folded-optics microscopes. This comparison will be included in Part 2.

Combined References and End Notes

- Beale, M.A. (1991) *The Lensman portable microscope--an aid to school fieldwork?* Journal of Biological Education (Society of Biology) Summer 1991, Vol 25, Issue 2
- Bracegirdle, Brian (1996) *Notes on Modern Microscope Manufacturers*. Oxford: Quekett Microscopical Club
- Dickinson, Rick (2013) Private emails. The author would like to thank Mr. Keith Dunning and Mr. Rick Dickinson for their valuable assistance in providing information and pictures. However, the content of this paper is the sole responsibility of the author
- Dingley, Mike (1998) *The Lensman Microscope: a review*. Micscape Magazine
<http://www.microscopy-uk.org.uk/mag/indexmag.html?http://www.microscopy-uk.org.uk/mag/art98/lensman.html>
- Dunning, Keith (See Dickinson, 2013)
- Enhelion (2002) *Micron by Enhelion*. Single sheet of instructions packaged with the microscope. July 2002
- Flickr (undated-1)
<http://www.flickr.com/photos/9574086@N02/8757351327/sizes/l/in/set-72157624402374235/>, and
<http://www.flickr.com/photos/9574086@N02/8866715919/in/set-72157624402374235/>
- Journal of Microscopy (1996) Published on behalf of the Royal Microscopical Society. Volume 183, Issue 2, pps 181-186, *Obituary for John Norrie McArthur 1901-1996*, Wiley-Blackwell, August 1996. Available on-line at
<http://onlinelibrary.wiley.com/doi/10.1046/j.1365-2818.1996.1050648.x/abstract>
- Kreindler, R. Jordan and Yuval Goren (2011-1) *Comparison of the Swift FM-31 Portable Field Microscope and an FM-31 Clone*. Micscape Magazine. March 2011
http://www.microscopy-uk.org.uk/mag/artmar11/A_Comparison_of_FM-31_and_clone.pdf

Kreindler, R. Jordan and Yuval Goren (2011-2) *The TWX-1 Folded-Optics Microscope: Monarch of the Folded-Optics Kingdom*. Micscape Magazine. November 2011

<http://www.microscopy-uk.org.uk/mag/artnov11/TWX-1.pdf>

Nikon (2013) Model H Microscope, MicroscopyU, The Source of Microscopy Education.

<http://www.microscopyu.com/museum/modelh.html>

Rowntree, Frank (1996)

<http://www.microscopy-uk.org.uk/mag/indexmag.html?http://www.microscopy-uk.org.uk/mag/articles/obituar1.html>

Saitoh, Motohiko (2006) McArthur microscope and Japan,

<http://www.geocities.jp/saitohmoto/essay/mcarthur/mcarthur-en.html>

The Herald, Glasgow (1996) *Obituary of Dr John Norrie McArthur, MRCS, LRCP, Hon FRMS, Hon Doctorate Open University, Hon FRS, born Glasgow, November 9, 1901, died Cambridge, April 26, 1996*. May 11, 1996

Walker, David (2005) *Microscopy on the move: A look at the Trekker field microscope*. Micscape Magazine. June 2005

©2013 Text and photographs (except as noted) by the author.

The author welcomes any suggestions for corrections or improvements. He has a continuing interest in folded-optics and field portable microscopes, and stereo and unique historical microscopes.

The author can be reached at:

R. Jordan Kreindler: leona111rjk@gmail.com

Published in the online magazine Micscape, July 2013

www.microscopy-uk.org.uk

Please report any Web problems or offer general comments to the Micscape Editor.

Micscape is the on-line monthly magazine of the Microscopy UK web
site at
Microscopy-UK